
mobi pro flexi









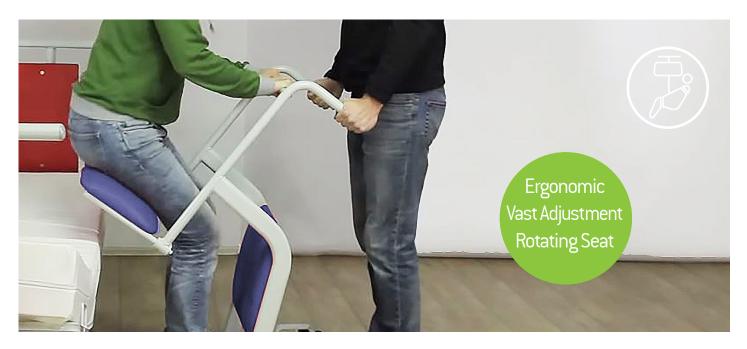

The Mobi Pro Flexi is an ergonomic and compact mobile assistive mobility device ideal for people with low leg strength and some trunk balance who experience difficulty walking. The Mobi Pro Flexi is an excellent safe patient handling solution for early mobilization. It offers support while standing and is an efficient device for patient transfers, improving safety and efficiency.

Rotating Seat

The Mobi Pro Flexi is an efficient and safe tool for the healthcare provider to move a patient, improving safety and caregiver time management. It features rotating padded seats and knee support for easy application. The feet are placed on a solid, firm footplate and then the patient pulls himself/herself up to a standing position using the integrated solid bracket. By simply rotating the seat cushion, the Mobi Pro becomes a comfortable seat that the patient can rest and be moved on.

Versatile Design

Using the Mobi Pro Flexi not only encourages mobilization and activity but gives the caregiver fast, stable and safe transport of a

patient to the bed, chair or the toilet and back if needed. The Mobi Pro Flexi is equipped with a mechanical mechanism to spread the legs, ensuring it will work with most seating (like a toilet or wheelchair). The vast adjustment range is easily set to three settings with a lever.

Quality

Our products are designed and manufactured to meet all the requirements of relevant standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B 117 43 Stockholm Phone:- +46 8 665 35 00 Fax:- +46 8 665 35 10 info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street Port Melbourne, VIC 3207 Phone: +61 1300 132 736 Fax: +61 [03] 9645 4954 info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road Ottawa, ON K2E 7K1 Phone: +613.723.6734 Fax: +613.723.1058 info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124 8076 PA Vierhouten Phone: +31 577 412 171 Fax: +31 577 412 170 info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K Austin, TX 78754 Phone: +512.476.7199 Fax: +512.476.7190 info.us@humancaregroup.com

Accessories



Foot Plate

Included with device. Durable ABS plastic. Replacement Part no. 1590035



Casters

Double wheel 7.5cm/3". Front Part no. 1690060 Rear with Brake Part no. 1690061



Ergonomic Grips

Made of rubber. Sold in pairs. Replacement Part no. 1510054

Lifting Situations



Use to/from a bed

Use to raise/lower a patient from/to a bed.



Use to/from a chair

Use to raise/lower a patient from/to a chair.



Use to/from the toilet

Use to raise/lower a patient from/to a toilet.



Tested to all relevent

standards, including ISO

14971, EN12182, MDD

93/42/EEG.

Features



Easy Movement

The Mobi Pro Flexi is easily moved by a single caregiver. A small turning radius makes it ideal for smaller spaces.



Adjustable Base

Manual adjustment of legs allows it to fit under most seating.



Padded Rotating Seats

Allows device to be used for seating as an alternative to a wheelchair.



Padded Knee Pad

Soft, cushioned knee pad to add patient comfort and reduce and device pressure at the knees while standing.

Technical Specifications

Mobi Pro Flexi			
Part No.	92500	Adjustable Leg (Min - Max)	43cm - 70cm / 17" - 27.5"
Width	58.5cm/23"	Handle Height (Inclined at Angle)	92cm - 100cm / 36.2" - 39.4"
Length	97cm/38.2"	Weight	19.6kg/43.2lb
Height	104cm/41"	Weight Capacity	150kg/330.6lb
Wheel Diameter	7.5cm/3"		

This product is manufactured in accordance with all relevant ISO standards.



Learn More About This Product



Follow Us on Social Media

