
niklas









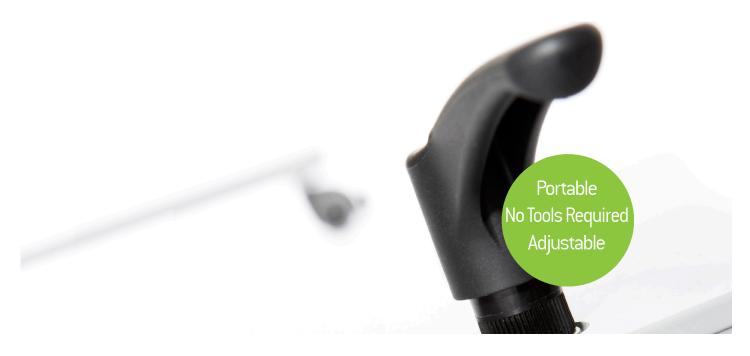

Our Niklas freestanding gantry system offers a unique way to safely and easily move patients without strain to the patient or caregiver. Designed for use with our series of ceiling lifts, the Niklas is portable and can be disassembled for easy transport. The Niklas is adjustable in height and width and requires no tools for assembly.

Portable Solution

The Niklas freestanding gantry system can easily be assembled without the use of tools. Unassembled, it can easily be transported in the carry case to be used between different locations like a home and cottage, or even in hotel rooms. It features dual wheel casters for maneuvrability and strong handles for easy movement. The adaptability of this gantry makes it an ideal option for temporary needs or rentals.

Adjustable

The unique design of the Niklas gantry allows it to be adjusted in both height and width, making it an extremely versatile gantry for moving patients from a bed, wheelchair, toilet or floor. The casters all lock into position for

added safety while the gantry is in use. The Niklas comes complete with a rail roller so it is ready for use with our portable lift systems.

Quality

Our products are designed and manufactured to meet all the requirements of relevant ISO, UL, IEC and EN standards. At every stage of the development process, every single product is given extensive attention in design, development and manufacturing. This means that a user of this product feels confident with its reliability and the caregiver feels confident with its safety.

HUMAN CARE SWEDEN (HQ)

Årstaängsvägen 21B 117 43 Stockholm Phone:- +46 8 665 35 00 Fax:- +46 8 665 35 10 info.se@humancaregroup.com

HUMAN CARE AUSTRALIA

Unit 4, 49 Bertie Street Port Melbourne, VIC 3207 Phone: +61 1300 132 736 Fax: +61 [03] 9645 4954 info.au@humancaregroup.com

HUMAN CARE CANADA

10-155 Colonnade Road Ottawa, ON K2E 7K1 Phone: +613.723.6734 Fax: +613.723.1058 info.ca@humancaregroup.com

HUMAN CARE NEDERLAND

Elspeterweg 124 8076 PA Vierhouten Phone: +31 577 412 171 Fax: +31 577 412 170 info.nl@humancaregroup.com

HUMAN CARE UNITED STATES

8006 Cameron Road, Suite K Austin, TX 78754 Phone: +512.476.7199 Fax: +512.476.7190 info.us@humancaregroup.com

Accessories



Rail Rollers

Optional rail rollers can be used with this product.
Niklas comes complete with a loop rail roller (shown).



Carry Case

Ideal for transport of gantry system.
Part no. 50599



Distance Strap

Reinforced nylon thread and metal hooks. Extend lift strap up to 100cm/39.3". Part no. 20124-xx (specify length)



Scale

Medical grade scale can be used on all Human Care lifts. Part no. 50590 (Roomer S), 50591 (Altair/ HeliQ), 50345 (NA).





Use to/from a bed

Can be placed over bed for transfers and repositioning with a ceiling lift/sling.



Use to/from a chair

Can be placed over a chair for transfers with a ceiling lift/sling.



Use to/from the floor

Can be placed over a patient in the supine position for transfers/lifting with a ceiling lift/sling.



Use to/from the toilet

Can be placed over a toilet for transfers with a ceiling lift/sling.

Features



Rail: Powder coated extruded aluminum Poles and feet: Powder coated steel.

Powder Coated



Portable Lifts

Compatible for use with our Altair series of ceiling lifts (shown) and our Roomer S.



Stationary Lifts

Compatible for use with our HeliQ series of ceiling lifts.



Tested to all relevent standards, including ISO 10535, hoists for the transfer of disabled persons.

Technical Specifications

Part Number	Description	Length	Height	Width	Product Weight	Weight Capacity
50600	Niklas	161.5cm to 262cm/ 63.6" to 103.2"	158cm to 238.8cm/ 62" to 94"	114.5cm / 45.1"	30kg/66lb	220kg/485lb

This product is manufactured in accordance with all relevant ISO standards.





Follow Us on Social Media



